

---

Product Manual


# 96-Well Cellular Senescence Assay Kit (SA- $\beta$ -gal Activity, Fluorometric Format)

## Catalog Number

CBA-231	120 assays
CBA-231-5	5 x 120 assays

**FOR RESEARCH USE ONLY**  
Not for use in diagnostic procedures

---

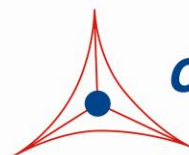

**CELL BIOLABS, INC.**  
*Creating Solutions for Life Science Research*

## **Introduction**

Normal primary cells proliferate in culture for a limited number of population doublings prior to undergoing terminal growth arrest and acquiring a senescent phenotype. This finite life span correlates with the age of the organism and with the life expectancy of the species from which the cells were obtained; such that the older the age or the shorter the life span, the less the ability of the cells to undergo population doubling. Senescent cells are characterized by an irreversible G<sub>1</sub> growth arrest involving the repression of genes that drive cell cycle progression and the upregulation of cell cycle inhibitors like p16<sup>INK4a</sup>, p53, and its transcriptional target, p21<sup>CIP1</sup>. They are resistant to mitogen-induced proliferation, and assume a characteristic enlarged, flattened morphology. Research into the pathways that positively regulate senescence and ways cells bypass senescence is therefore critical in understanding carcinogenesis. Normal cells have several mechanisms in place to protect against uncontrolled proliferation and tumorigenesis.

Senescent cells show common biochemical markers such as expression of an acidic senescence-associated  $\beta$ -galactosidase (SA- $\beta$ -Gal) activity. While senescence has been characterized primarily in cultured cells, there is also evidence that it occurs *in vivo*. Cells expressing markers of senescence such as SA- $\beta$ -Gal have been identified in normal tissues.

The 96-well Cellular Senescence Assay Kit provides an easy-to-use and efficient method to determine the cellular senescence by measuring SA- $\beta$ -Gal activity using a fluorometric substrate. This quantitative assay uses cell lysate for both SA- $\beta$ -galactosidase activity determination and normalization of samples containing different cell numbers. Each kit provides sufficient quantities to perform up to 120 assays in a 96-well plate.

## **Related Products**

1. CBA-230: Cellular Senescence Assay Kit (SA- $\beta$ -gal Staining)
2. CBA-232: Quantitative Cellular Senescence Assay (SA  $\beta$ -Gal)
3. CBA-240: CytoSelect™ Cell Viability and Cytotoxicity Assay
4. AKR-100:  $\beta$  Galactosidase Staining Kit

## **Kit Components (shipped at room temperature)**

1. 2X Cell Lysis Buffer (Part No. 123101): One 10 mL bottle
2. 2X Reaction Buffer (Part No. 123102): One 10 mL bottle
3. SA- $\beta$ -Gal Substrate (20X) (Part No. 123103): One 300  $\mu$ L amber tube
4. Stop Solution (Part No. 123104): One 25 mL bottle

## **Materials Not Supplied**

1. Senescent cells or tissue samples
2.  $\beta$ -mercaptoethanol
3. 96-well plate suitable for a fluorescence plate reader
4. Protein Assay Reagents

## **Storage**

Store SA- $\beta$ -gal substrate solution protected from light at -20°C. Store all other components at room temperature.

## **Preparation of Reagents**

- 1X Cell Lysis Buffer: Prepare a 1X Cell Lysis Buffer by diluting the provided 2X stock 1:2 in ddH<sub>2</sub>O. Store the diluted solution at room temperature for up to six months. Immediately before use, add proper amount of proteinase inhibitors such as PMSF.
- 2X Assay Buffer: Immediately before use, add  $\beta$ -mercaptoethanol to 2X Reaction Buffer at a final concentration of 10 mM and dilute 20X SA- $\beta$ -Gal Substrate to 1X with 2X Reaction Buffer containing 10 mM  $\beta$ -mercaptoethanol. Don't store 2X Assay Buffer.


<b>Reagent</b>	<b>96-well</b>	<b>24-well</b>	<b>6-well</b>	<b>10 cm</b>
1X Cell Lysis Buffer	100 $\mu$ L	400 $\mu$ L	1000 $\mu$ L	1500 $\mu$ L

## **Assay Protocol**

1. Aspirate the medium from the senescent cells.
2. Wash the cells once with 200  $\mu$ L of cold 1X PBS and aspirate the wash.
3. Add 100  $\mu$ L of cold 1X Cell Lysis Buffer (see the table above for the required amount of 1X Cell Lysis Buffer of other plate formats). Incubate at 4°C for 5 minutes. Transfer the whole lysate to a microcentrifuge tube and centrifuge 10 minutes at 4°C. Collect supernatant as cell lysate.
4. (optional) Determine the total protein concentration of each cell lysate sample by protein assay such as Pierce's BCA protein Assay.
5. Transfer 50  $\mu$ L of the cell lysate to a 96-well plate. Add 50  $\mu$ L of freshly prepared 2X Assay Buffer. Incubate the wells at 37°C protected from light for 1- 3 hr.
6. Remove 50  $\mu$ L of the reaction mixture to a 96-well plate suitable for fluorescence measurement. Stop the reaction by adding 200  $\mu$ L of Stop solution.
7. Read fluorescence with a fluorescence plate reader at 360 nm (Excitation) / 465 nm (Emission).

## **Example of Results**

The following figures demonstrate typical with the 96-well Cellular Senescence Assay Kit. One should use the data below for reference only. This data should not be used to interpret actual results.

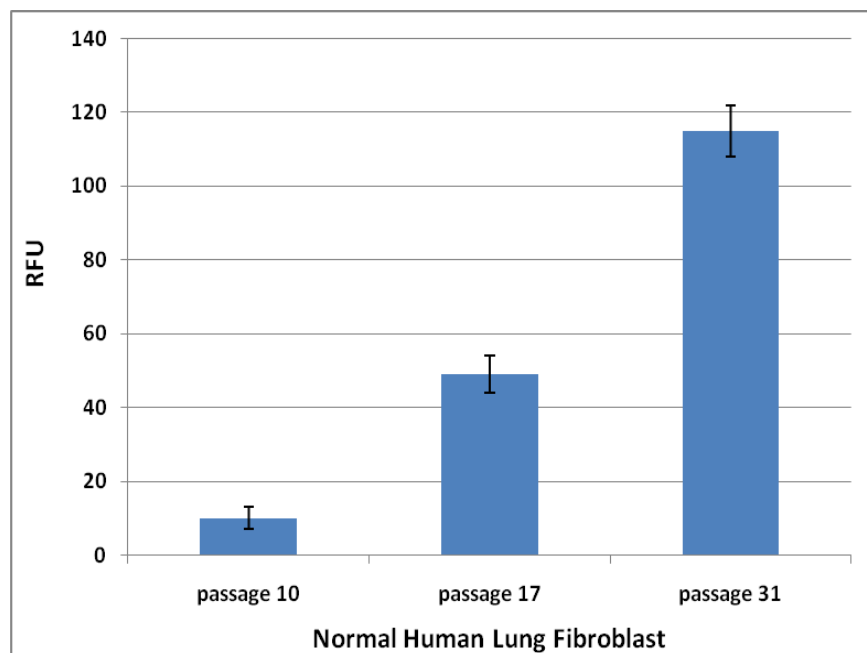

**Figure 1: SA-β-Gal activity in Senescent Human Lung Fibroblast HFL-1 Cells.** Normal HFL-1 cells with different passage numbers were lysed. Lysates were allowed to incubate with SA-β-Gal Substrate for 1 hr at 37°C. SA-β-Gal activities were measured as described in the Assay Protocol.

## References

1. Current Protocols in Molecular Biology, John Wiley & Sons Press.
2. Campisi, J. (2000) *In Vivo* 14, 183-188.
3. Dimri, G. P., X. Lee, G. Basile, M. Acosta, G. Scott, C. Roskelley, E. E. Medrano, M. Linskens, I. Rubelj, O. Pereira-Smith, M. Peacocke, and J. Campisi. (1995) *Proc. Natl. Acad. Sci. USA* 92:9363-9367.

## Recent Product Citations

1. Li, Z. et al. (2026). Gastrodia elata-Derived Parishin Extracts Attenuate Aging by Modulating Oxidative Stress, Inflammation, Apoptosis, and Metabolism. *Food Sci Nutr.* 14(4):e71737. doi: 10.1002/fsn3.71737.
2. Tsutsui, Y. et al. (2026). Human iPSC-based Modeling of Pulmonary Fibrosis Reveals p300/CBP Inhibition Suppresses Alveolar Transitional Cell State. *Nat Commun.* 17(1):1214. doi: 10.1038/s41467-026-68909-z.
3. Zhou, Y. et al. (2025). Antioxidant and Anti-Senescence Polyvinyl Alcohol-Gallic Acid Supramolecular Hydrogels for Stem Cell Culture. *Adv Healthc Mater.* doi: 10.1002/adhm.202402882.
4. Tarasova, K. et al. (2025). Dexamethasone: a double-edged sword in the treatment of osteoarthritis. *Sci Rep.* 15(1):11832. doi: 10.1038/s41598-025-96050-2.
5. Dzubanov, M. et al. (2025). Protocol for isolation of human bone marrow stromal cells and characterization of cellular metabolism. *STAR Protoc.* 6(1):103553. doi: 10.1016/j.xpro.2024.103553.

6. Haroon, M. & Kang, S.C. (2024). Kaempferol Synergistically Enhances Cisplatin-induced Apoptosis and Cell Cycle Arrest in Colon Cancer Cells. *J Cancer Prev.* **29**(3):69-87. doi: 10.15430/JCP.24.013.
7. Namirah, I. et al. (2024). The effect of ethanol-based coriander (*Coriandrum sativum* L.) seed extract on oxidative stress, antioxidant level and cellular senescence in the heart of obese rat. *J Pharm Pharmacogn Res.* **12**(6):1111-1120. doi: 10.56499/jppres24.1957\_12.6.1111.
8. Zhang, L. et al. (2024). Nesfatin-1 attenuated lipopolysaccharide-induced inflammatory response and senescence in human dental pulp cells. *Heliyon.* **10**(12):e32108. doi: 10.1016/j.heliyon.2024.e32108.
9. Boulestreau, J. et al. (2024). Extracellular vesicles from senescent mesenchymal stromal cells are defective and cannot prevent osteoarthritis. *J Nanobiotechnology.* **22**(1):255. doi: 10.1186/s12951-024-02509-1.
10. Oram, M.K. et al. (2024). RNF4 prevents genomic instability caused by chronic DNA under-replication. *DNA Repair (Amst).* doi: 10.1016/j.dnarep.2024.103646.
11. Hardiany, N.S. et al. (2024). Exploration of neuroprotective effect from *Coriandrum sativum* L. ethanolic seeds extracts on brain of obese rats. *Sci Rep.* **14**(1):603. doi: 10.1038/s41598-024-51221-5.
12. Bai, W. et al. (2024). S-nitrosylation of AMPK $\gamma$  impairs coronary collateral circulation and disrupts VSMC reprogramming. *EMBO Rep.* **25**(1):128-143. doi: 10.1038/s44319-023-00015-3.
13. Benova, A. et al. (2023). Omega-3 PUFAs prevent bone impairment and bone marrow adiposity in mouse model of obesity. *Commun Biol.* **6**(1):1043. doi: 10.1038/s42003-023-05407-8.
14. Shen, C.Y. et al. (2023). Unveiling the molecular basis of inflamm-aging induced by advanced glycation end products (AGEs)-modified human serum albumin (AGE-HSA) in patients with different immune-mediated diseases. *Clin Immunol.* **252**:109655. doi: 10.1016/j.clim.2023.109655.
15. Ambrosio, M.R. et al. (2023). Targeting G-quadruplex motifs interferes with differentiation of adipose-derived mesenchymal stem cells. *Stem Cell Res Ther.* **14**(1):98. doi: 10.1186/s13287-023-03320-9.
16. Choi, D. et al. (2023). Vertical Vibration of Mouse Osteoblasts Promotes Cellular Differentiation and Cell Cycle Progression and Induces Aging In Vitro. *Biomedicines.* **11**(2):444. doi: 10.3390/biomedicines11020444.
17. Takami, Y. et al. (2022). Novel pathophysiological roles of  $\alpha$ -synuclein in age-related vascular endothelial dysfunction. *FASEB J.* **36**(10):e22555. doi: 10.1096/fj.202101621R.
18. Nelczyk, A.T. et al. (2022). The nuclear receptor TLX (NR2E1) inhibits growth and progression of triple-negative breast cancer. *Biochim Biophys Acta Mol Basis Dis.* **1868**(11):166515. doi: 10.1016/j.bbadis.2022.166515.
19. Chelyapov, N. et al. (2022). Autologous NK cells propagated and activated ex vivo decrease senescence markers in human PBMCs. *Biochem Biophys Rep.* doi: 10.1016/j.bbrep.2022.101380.
20. Furukawa, M. et al. (2022). Analysis of senescence in gingival tissues and gingival fibroblast cultures. *Clin Exp Dent Res.* doi: 10.1002/cre2.581.
21. Kelleher, A.M. et al. (2021). Deficiency of PARP-1 and PARP-2 in the mouse uterus results in decidualization failure and pregnancy loss. *Proc Natl Acad Sci U S A.* **118**(40):e2109252118. doi: 10.1073/pnas.2109252118.
22. Deng, Z. et al. (2021). Biofunction of Polydopamine Coating in Stem Cell Culture. *ACS Appl Mater Interfaces.* doi: 10.1021/acsami.0c22565.
23. Lee, Y.N. et al. (2021). Ultrasonic microbubble VEGF gene delivery improves angiogenesis of senescent endothelial progenitor cells. *Sci Rep.* **11**(1):13449. doi: 10.1038/s41598-021-92754-3.

24. Hirata, Y. et al. (2021). Advanced maternal age induces fetal growth restriction through decreased placental inflammatory cytokine expression and immune cell accumulation in mice. *J Reprod Dev.* doi: 10.1262/jrd.2021-034.
25. Mehdi, S.J. et al. (2021). Normal and cancer fibroblasts differentially regulate TWIST1, TOX and cytokine gene expression in cutaneous T-cell lymphoma. *BMC Cancer.* **21**(1):492. doi: 10.1186/s12885-021-08142-7.
26. Bourdon, B. et al. (2021). Marine Collagen Hydrolysates Promote Collagen Synthesis, Viability and Proliferation While Downregulating the Synthesis of Pro-Catabolic Markers in Human Articular Chondrocytes. *Int. J. Mol. Sci.* **22**(7):3693. doi: 10.3390/ijms22073693.
27. Kong, C.S. et al. (2021). Embryo biosensing by uterine natural killer cells determines endometrial fate decisions at implantation. *FASEB J.* **35**(4):e21336. doi: 10.1096/fj.202002217R.
28. Baxley, R.M. et al. (2021). Bi-allelic MCM10 variants associated with immune dysfunction and cardiomyopathy cause telomere shortening. *Nat Commun.* **12**(1):1626. doi: 10.1038/s41467-021-21878-x.
29. Yamaguchi, S. et al. (2021). Characterization of an active LINE-1 in the naked mole-rat genome. *Sci Rep.* **11**(1):5725. doi: 10.1038/s41598-021-84962-8.
30. Bourdon, B. et al. (2021). Marine Collagen Hydrolysates Downregulate the Synthesis of Pro-Catabolic and Pro-Inflammatory Markers of Osteoarthritis and Favor Collagen Production and Metabolic Activity in Equine Articular Chondrocyte Organoids. *Int J Mol Sci.* **22**(2):E580. doi: 10.3390/ijms22020580.

## **Warranty**

These products are warranted to perform as described in their labeling and in Cell Biolabs literature when used in accordance with their instructions. THERE ARE NO WARRANTIES THAT EXTEND BEYOND THIS EXPRESSED WARRANTY AND CELL BIOLABS DISCLAIMS ANY IMPLIED WARRANTY OF MERCHANTABILITY OR WARRANTY OF FITNESS FOR PARTICULAR PURPOSE. CELL BIOLABS' sole obligation and purchaser's exclusive remedy for breach of this warranty shall be, at the option of CELL BIOLABS, to repair or replace the products. In no event shall CELL BIOLABS be liable for any proximate, incidental or consequential damages in connection with the products.

## **Contact Information**

Cell Biolabs, Inc.  
5628 Copley Drive  
San Diego, CA 92111  
Worldwide: +1 858-271-6500  
USA Toll-Free: 1-888-CBL-0505  
E-mail: [tech@cellbiolabs.com](mailto:tech@cellbiolabs.com)  
[www.cellbiolabs.com](http://www.cellbiolabs.com)

©2007-2026: Cell Biolabs, Inc. - All rights reserved. No part of these works may be reproduced in any form without permissions in writing.