

Buffers and Stabilizers



PROVIDING CONSCIOUS IMMUNOASSAY SOLUTIONS

kementec

Distributed by
cliniSciences


Buffers and Stabilizers

Kementec Solutions' ECO-TEK buffers are designed for use in colorimetric as well as in chemiluminescence based assays such as ELISA, Western Blotting and IHC. They are all ready-to-use and robust diluents offering superior stabilization of proteins conjugated to Horseradish Peroxidase (HRP) and Alkaline Phosphatase (AP). These stabilizers prevent loss of enzyme activity and maintain the proteins' molecular conformation. Furthermore, all our ECO-TEK buffers are Bovine Serum Albumin (BSA) free and several are animal-free.

WellChampion®

Plate Blocking & Protection Buffer for pre-coated plates

WellChampion preserves and protects biological activity of proteins and other biomolecules coated and dried on microwell surfaces. The protection of the coated and dried surface prevents degradation, denaturation and leaching. As WellChampion also efficiently blocks any free binding sites left on the microwell surface after coating, separate blocking of the surface is not necessary. WellChampion dramatically speeds up your plate coating and stabilization process in most cases at the same time resulting in improved assay sensitivity.

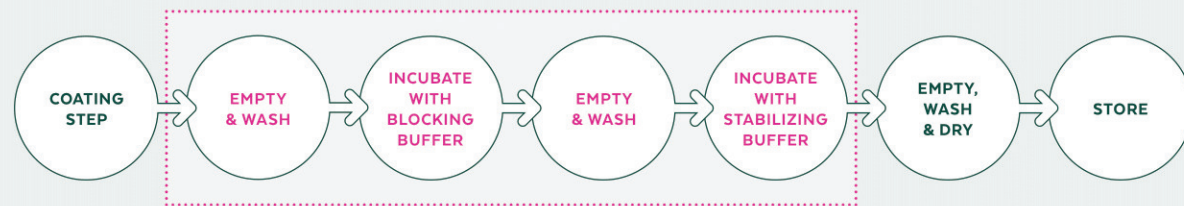

- Coat the microwell plate and add WellChampion directly, without emptying the plate first.
- Allow the surface to saturate with WellChampion.
- Empty the plate without tapping and let it dry.
- The plate is now blocked, providing long-term stability.



Choose the right buffer

The following table will guide you to the optimal choice of our different solutions.

D: Diluent, S: Stabilizer, B: Blocker

Animal-free Buffers	Standards & Controls	Samples	HRP-conjugates	AP-conjugates	Antibodies	Surfaces
HRP-StabilPLUS *			D/S			
AP-StabilPLUS				D/S		
UNI-StabilPLUS			D/S	D/S		
Protein-StabilPLUS	D/S				D/S	
Antibody Enhancer *			D/S	D/S	D/S	
Sample X-tra	D	D				
Effect Diluent		D	D	D	D	
Synthetic Blocker						B
WellChampion® *						B/S

* Not Animal-free

Synthetic Blocking Buffers



Blocking Buffers for ELISA and Blotting

Purely synthetic, protein free blocking buffers for ELISA and Blotting applications with high blocking capacity: Improving signal-to-noise ratio and eliminating the risk of false positive cross-reactions in the assay.

HRP-StabilPLUS

HRP Conjugate Diluent & Stabilizer

HRP-StabilPLUS is our most successful and widely used conjugate buffer for HRP labeled proteins due to its unsurpassed long-lasting stabilizing effect in both high and low protein concentrations.

AP-StabilPLUS



AP Conjugate Diluent & Stabilizer

AP-StabilPLUS is a highly recommended alternative to conventional Alkaline Phosphatase conjugate stabilizers. AP-StabilPLUS provides a long-lasting stabilizing effect.

UNI-StabilPLUS



AP & HRP Conjugate Diluent & Stabilizer

UNI-StabilPLUS is a universal conjugate diluent designed for the stabilization of both AP- & HRP-labelled proteins. Utilizing the same buffer for all conjugates simplifies daily work processes.

Antibody Enhancer

HRP- & AP- Antibody Diluent, Stabilizer & Amplifier

Antibody Enhancer is an antibody efficient conjugate stabilizer as well as signal amplifier. The buffer enhances the binding capacity of the diluted antibodies and provides the highest assay signal-to-noise ratio among our stabilizers.

Sample X-tra



Sample dilution buffer for Immunoassays

The dilution of sample matrices & assay standards with Sample X-tra minimizes nonspecific interactions between the matrix proteins and the surface of the plate, improving signal-to-noise ratio.

Effect Diluent buffers, Animal-free



Sample dilution buffers for Immunoassays

The 3 Effect Diluent buffers (Low, Medium and High) eliminate interference from HAMAs (Human-Anti-Mouse Antibodies), Heterophilic Abs and RFs (Rheumatoid Factors) by blocking cross-reactions without affecting the correct binding of the assay antibodies.

Protein-StabilPLUS



Protein Stabilizer

Protein-StabilPLUS gives superior stabilization of ready-to-use pre-diluted standards and controls used in immunoassays.



Animal-free



Our immunoassay solutions are developed in order to create a healthy and safe work environment and preserve natural resources, while helping our customers fulfill significant regulatory requirements.

We call it ECO-TEK®

Profile

Kementec is a well-established, innovative, and service oriented company with more than 35 years of experience as a manufacturer and supplier of components to the industry and to researchers worldwide.


Being a market leader in developing solutions for immunoassays, our continuous goal is to deliver sustainable products without compromising on quality, sensitivity and accuracy.

Kementec values long lasting business relationships and the success of our customers. We consider ourselves more as a trusted partner than as a general supplier. A part of this vision is providing high quality colorimetric substrates and color-coded reagents, described as:

Colors of Success!

How to navigate through our website?

Homepage: When you first visit our website, you'll land on the homepage. This is where you can find an overview of our services, products, or any important announcements. From here, you can navigate to other sections of the website.

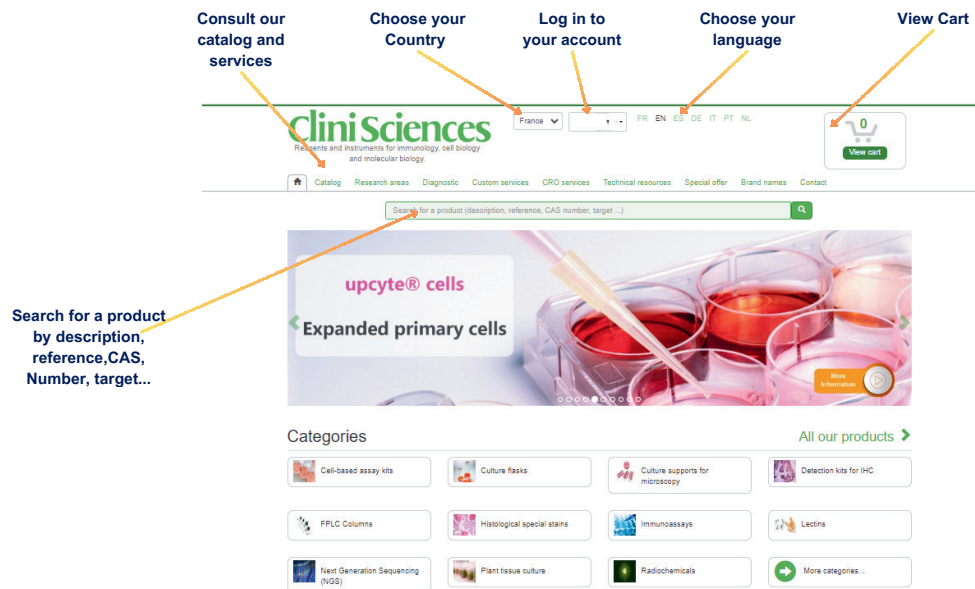

Search Bar: Look for specific Products and equipment

Search by selecting either RUO or IVD products

Refine your search by choosing from the following options

Apply Filter

Search for a product (description, reference, CAS number, target ...)

Search result : 39126007 product found

Cat#	Description	Cond.
NB-45-00042-100	Super Ni-NTA Agarose Resin	100ml
NB-45-00042-25	Super Ni-NTA Agarose Resin	25ml
NB-45-00058-4	Proteus 1-step Batch Mini Spin Column Pack	40pc
NB-12-6001-3	NeoLine pipette 2-20 µl	1unit
NB-12-0023C	Mini Centrifuge N500C @10,000rpm (including 6x1.5/2.0ml angle rotor)	1pcs
NB-03-0160	Proteinase K (Powder)	100mg
NB-60-0001	NeoPrep mini	50columns
NB-12-8001-19	Combs for NeoPRO4 mini (1.5mm, 15 wells)	5pieces
NB-12-8001-20	Spacer glasses flat for NeoPRO4 mini (0.75mm, 100*73mm)	5pieces
NB-12-8001-04	Short glasses flat for NeoPRO4 mini (1.0mm, 100*73mm)	10pieces

Refine your search:

☒ RUO ☒ CE / IVD

Conjugate (100) +

Reactivity species (100) +

Class (66) +

Host species (100) +

Cell name (2) +

Application (100) +

Clone (100) +

APPLY FILTERS

REINITIALIZE

SHARE THIS SEARCH

“Remember, if you ever get lost or have questions, don’t hesitate to reach out to our customer support team. They’re here to help you address any concerns you may have.”

CliniSciences

Austria

Company: CliniSciences GmbH
Address: Sternwartestrasse 76, A-1180
Wien - Austria
Telephone: +43 720 115 580
Fax: +43 720 115 577
Email: österreich@clinisciences.com
Web: <https://www.clinisciences.com>



Belgium

Company: CliniSciences S.R.L
Address: Avenue Stalingrad 52, 1000
Brussels - Belgium
Telephone: +32 2 31 50 800
Fax: +32 2 31 50 801
Email: belgium@clinisciences.com
Web: <https://www.clinisciences.com>



Denmark

Company: CliniSciences ApS
Address: Oesterbrogade 226, st. 1,
Copenhagen, 2100 - Denmark
Telephone: +45 89 888 349
Fax: +45 89 884 064
Email: denmark@clinisciences.com
Web: <https://www.clinisciences.com>



Finland

Company: CliniSciences ApS
Address: Oesterbrogade 226, st. 1,
Copenhagen, 2100 - Denmark
Telephone: +45 89 888 349
Fax: +45 89 884 064
Email: suomi@clinisciences.com
Web: <https://www.clinisciences.com>



France

Company: CliniSciences S.A.S
Address: 74 Rue des Suisses, 92000
Nanterre- France
Telephone: +33 9 77 40 09 09
Fax: +33 9 77 40 10 11
Email: info@clinisciences.com
Web: <https://www.clinisciences.com>



Germany

Company: Clinisiences GmbH
Address: Wilhelm-Mauser-Str. 41-43,
50827 Köln - Germany
Telephone: +49 221 9498 320
Fax: +49 221 9498 325
Email: deutschland@clinisciences.com
Web: <https://www.clinisciences.com>



Iceland

Company: CliniSciences ApS
Address: Oesterbrogade 226, st. 1,
Copenhagen, 2100 - Denmark
Telephone: +45 89 888 349
Fax: +45 89 884 064
Email: island@clinisciences.com
Web: <https://www.clinisciences.com>



Ireland

Company: CliniSciences Limited
Address: Ground Floor, 71 lower Baggot street
Dublin D02 P593 - Ireland
Telephone: +353 1 6971 146
Fax: +353 1 6971 147
Email: ireland@clinisciences.com
Web: <https://www.clinisciences.com>



Italy

Company: CliniSciences S.r.l
Address: Via Maremmana inferiore 378
Roma 00012 Guidonia Montecelio - Italy
Telephone: +39 06 94 80 56 71
Fax: +39 06 94 80 00 21
Email: italia@clinisciences.com
Web: <https://www.clinisciences.com>



Netherlands

Company: CliniSciences B.V.
Address: Kraijenhoffstraat 137A,
1018RG Amsterdam, - Netherlands
Telephone: +31 85 2082 351
Fax: +31 85 2082 353
Email: nederland@clinisciences.com
Web: <https://www.clinisciences.com>



Norway

Company: CliniSciences AS
Address: c/o MerVerdi Munkerudtunet 10
1164 Oslo - Norway
Telephone: +47 21 988 882
Email: norge@clinisciences.com
Web: <https://www.clinisciences.com>



Poland

Company: CliniSciences sp.Z.o.o.
Address: ul. Rotmistrza Witolda Pileckiego 67
lok. 200 - 02-781 Warszawa -Poland
Telephone: +48 22 307 0535
Fax: +48 22 307 0532
Email: polska@clinisciences.com
Web: <https://www.clinisciences.com>



Portugal

Company: Quimigen Unipessoal LDA
Address: Rua Almada Negreiros, Lote 5, Loja 14,
2615-275 Alverca Do Ribatejo - Portugal
Telephone: +351 30 8808 050
Fax: +351 30 8808 052
Email: info@quimigen.com
Web: <https://www.quimigen.pt>



Spain

Company: CliniSciences Lab Solutions
Address: C/ Hermanos del Moral 13
(Bajo E), 28019, Madrid - Spain
Telephone: +34 916 750 700
Fax: +34 91 269 40 74
Email: espana@clinisciences.com
Web: <https://www.clinisciences.com>



Sweden

Company: CliniSciences ApS
Address: Oesterbrogade 226, st. 1,
Copenhagen, 2100 - Denmark
Telephone: +45 89 888 349
Fax: +45 89 884 064
Email: sverige@clinisciences.com
Web: <https://www.clinisciences.com>



Switzerland

Company: CliniSciences AG
Address: Fracht Ost Flughafen Kloten
CH-8058 Zürich - Switzerland
Telephone: +41 (044) 805 76 81
Fax: +41 (044) 805 76 75
Email: switzerland@clinisciences.com
Web: <https://www.clinisciences.com>



UK

Company: CliniSciences Limited
Address: 11 Progress Business center, Whittle
Parkway, SL1 6DQ Slough- United Kingdom
Telephone: +44 (0)1753 866 511
or +44 (0) 330 684 0982
Fax: +44 (0)1753 208 899
Email: uk@clinisciences.com
Web: <https://www.clinisciences.com>



USA

Company: CliniSciences LLC
Address: c/o Carr Riggs Ingram,
500 Grand Boulevard, Suite 210 Miramar
Beach, FL 32550- USA
Telephone: +1 850 650 7790
Fax: +1 850 650 4383
Email: usa@clinisciences.com
Web: <https://www.clinisciences.com>

